

3M™ Tegaderm™ Transparent Film Dressing Frame Style

1626, 1627, 1628, 1629, 1630, 1634, 1622W, 1624W, 1626W, 9505W, 9506W

3M™ Tegaderm™ HP Transparent Film Dressing Frame Style

9519HP, 9534HP, 9536HP, 9543HP, 9545HP, 9546HP, 9548HP

3M™ Tegaderm™ Transparent Film Dressing with Border

1614, 1616

3M™ Tegaderm™ I.V. Transparent Film Dressing with Border

1610, 1650, 1655, 9525HP

3M™ Tegaderm™ Transparent Film Dressing First Aid Style

1620, 1621

GB

Product Description

Tegaderm™ Film consists of a thin film backing with a non-latex, hypoallergenic adhesive. Tegaderm™ Film with Border is notched and reinforced with soft cloth tape to provide a better seal around catheters and other devices. The dressing is breathable, allowing good oxygen and moisture vapor exchange. It is waterproof and impermeable to liquids, bacteria, and viruses.* An intact dressing protects the site from outside contamination.

*In vitro testing shows that the film of Tegaderm™ and Tegaderm™ HP dressings provides a viral barrier from viruses 27 nm in diameter or larger while the dressing remains intact without leakage.

Indications for Use

Tegaderm™ Film can be used to cover and protect catheter sites and wounds, to maintain a moist environment for wound healing or to facilitate autolytic debridement, as a secondary dressing, as a protective cover over at-risk skin, to secure devices to the skin, to cover first and second degree burns, and as a protective eye covering. Do not use the dressing as a replacement for sutures and other primary wound closure methods.

Warning

Using the Tegaderm dressing to cover a drug delivery patch could lead to faster absorption, increasing the risk for overdose.


Precautions

1. Stop any bleeding at the site before applying the dressing.
2. Do not stretch the dressing during application as tension can cause skin trauma.
3. Make sure the skin is clean, free of soap residue and lotion and allowed to dry thoroughly before applying the dressing to prevent skin irritation and to ensure good adhesion.
4. The dressing may be used on an infected site, only when under the care of a health care professional.
5. Antimicrobial ointments containing polyethylene glycols may compromise the strength of the Tegaderm™ HP Transparent Film Dressings.


Shelf Life: For shelf life, refer to the expiration date on the package.

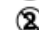
For further information outside of the United States, contact your local 3M representative or contact us at 3M.com and select your country.

Explanation of Symbols

 Not Made With Natural Rubber Latex

 Caution, see instructions for use


 Do not use if package is damaged

 Do not reuse

 Use by date

 Batch code

 Manufacturer

 Date of manufacture

 Sterilized using irradiation

CE 2797

Made in U.S.A. for
3M Health Care
2510 Conway Ave.
St. Paul, MN 55144 U.S.A.
1-800-228-3957
3M.com/Tegaderm

3M and Tegaderm are trademarks of 3M.
Used under license in Canada.
© 2020, 3M. All rights reserved.
Issue Date: 2020-08
34-8726-0592-7

3M Canada Company
P.O. Box/C.P. 5757
London, Ontario, N6A 4T1
1-800-364-3577
3M.com/Tegaderm

3M et Tegaderm sont des
marques de commerce de 3M,
utilisées sous licence au Canada.
© 2020, 3M. Tous droits réservés.

EC REP 3M Deutschland GmbH
Health Care Business
Carl-Schurz-Str. 1
41453 Neuss
Germany